
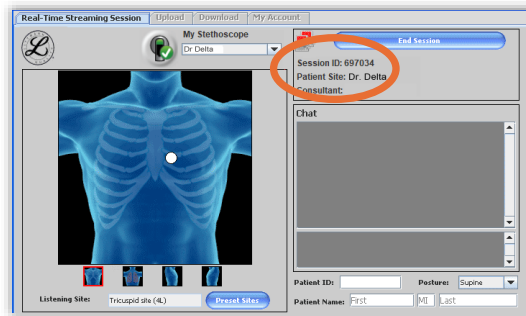
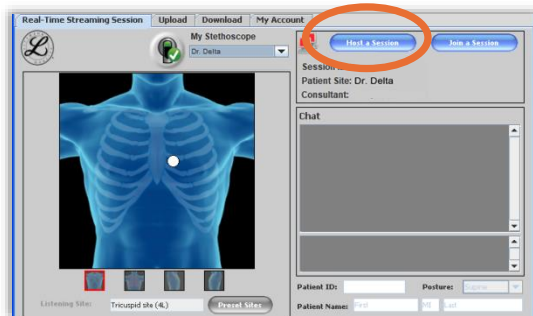




Quick Reference Guide – Patient site

1. To securely connect to the internet, start the **AnyConnect**  app and log in.
2. Turn on the **Telesteth Littman 3200 stethoscope** by pressing the Power button.
3. Start the **Telesteth** app.
 - A Java popup might appear briefly as the app starts up.
4. Enter your OTN Telesteth username and password.
5. Follow the instructions to connect your stethoscope with your computer.
 - Click on the My Stethoscope down arrow, pick your stethoscope from the list and follow the instructions in the pop-up window.
 - You will see a green checkmark when connected.
6. As the Patient site, click 'Host a Session'. A random Session ID number will be generated. Share the Session ID number with the Consultant site.




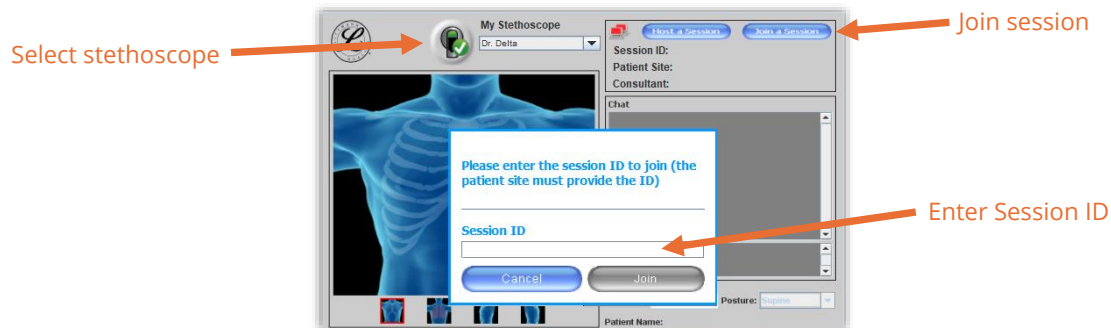
7. When successfully connected, the computer icon turns green . Follow the direction of the consultant regarding requested placement.




8. When the patient assessment is complete, click End Session and log out of the app.
9. Disconnect by right clicking the AnyConnect icon  in the taskbar and select Disconnect from the pop-up menu.


Quick Reference Guide – Consultant site

1. To securely connect to the internet, start the **AnyConnect**  app and log in.
2. Turn on the **Telesteth Littman 3200 stethoscope** by pressing the Power button.
3. Start the **Telesteth** app.
 - A Java popup might appear briefly as the app starts up.
4. Enter your OTN Telesteth username and password.
5. Follow the instructions to connect your stethoscope with your computer.
 - Click on the My Stethoscope down arrow, pick your stethoscope from the list and follow the instructions in the pop-up window.
 - You will see a green checkmark when connected.
6. As the Consultant site, click 'Join a Session' and enter the Session ID provided by the Patient site. (The Patient site will always 'Host a Session'.)



7. When successfully connected, the computer icon turns green  and you are ready to provide placement direction to the patient site. Adjust the filter settings and amplification using the buttons on the stethoscope.



8. When the patient assessment is complete, click End Session and log out of the app.
9. Disconnect by right clicking the AnyConnect icon  in the taskbar and select Disconnect from the pop-up menu.